

Boston
Scientific



Quick Start Guide

Charging System

SPECTRA
WAVEWRITER™

Spinal Cord Stimulator System



Ease of use is at the core of your SCS charging system

*The Charging System
aims to make life
as simple as
possible for you.*

An overview

Charger

Indicator Light



Power

+



Charging your device is as easy as 1-2-3!

- The charging system is **simple and easy to use**.
- The charging belt is **completely cordless**, so you can **enjoy full mobility** while charging your device.
- **You can decide** when and where you charge your device.



All you need to charge your rechargeable Spectra WaveWriter™ stimulator is:



A **power supply**



The **base station**



The **charger**




The **charging belt**



Charge your stimulator in 3 easy steps

1. Charge your charger

- Plug the base station into a power outlet and place the charger in the base station. A **yellow light** will appear on the charger. This means the **charging is in process**.

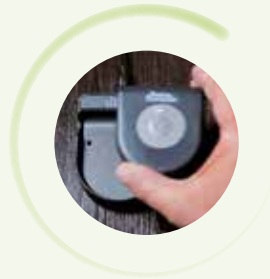

- Please note that as long as the **yellow light** shows, the charger is only **partially charged**. This means it can still be used to charge your stimulator but it may not be possible to recharge the stimulator fully.



- When a **green light** shows, the charger is **fully charged**. This means it is ready to charge your stimulator.




- If the indicator light is **not showing** while the charger is in the base station, there may be a **loose connection**. Check that the charger is securely placed in the base station.




2. Place the charger within your charging belt

- Once sufficiently charged, remove the charger from the base station. Before you place your charger in your charging belt, make sure to **switch your charger on**, using the circular **power button** in the middle of the charger.
- The **green indicator light** will show and the charger will **beep in a steady continuous rhythm**. Once the charger is correctly aligned with your stimulator, the **beeping will stop**.

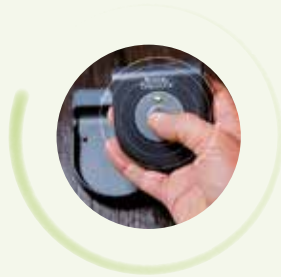


-
- Place the **charger** into the pouch of your charger belt and place the belt around your waist.
 - Adjust the straps on the belt so they are loose enough for you to feel comfortable. Make sure the **mesh fabric of your belt is facing outwards**.

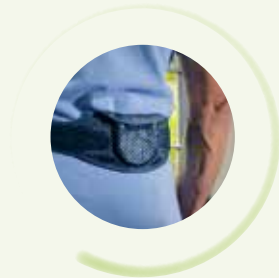

3. Position the belt to align the charger with your stimulator

- Move the belt slightly up and down and from side to side until the charger is **resting over your stimulator**. If you need to, you can adjust the length of the belt to the right size for you. When correctly positioned the **continuous single beeps will stop**.
- Check the charger occasionally to ensure it's **centered over your stimulator**. If not centered properly, charging times may vary.
- You know the stimulator is fully charged when you hear a **double beep for a total of 60 seconds** and see that the **green indicator light** has switched off. Although you don't have to charge your stimulator to 100% it can be helpful to know that your battery is fully charged.

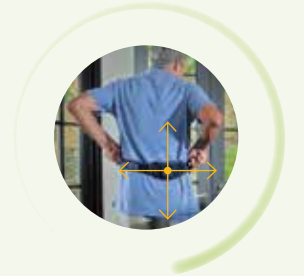

Tip: Once your stimulator is charged, return your charger to the base station. This ensures it will be fully charged for the next time you need it. Then you're all set!

Making charging work for you

Top Tips

- **Choose a charging schedule that works for you:** You can decide whether to have shorter, more frequent charging sessions or longer, less frequent charging sessions. Pick a routine that best fits into your daily lifestyle and don't forget that you can easily charge while watching your favorite TV program or reading the morning paper. A standard charging session once a week generally takes **30–60 minutes**.
- **Check your battery status via your Freelink™ Remote Control:** As long as your charger isn't currently charging your stimulator, your Freelink™ remote will show your stimulator's battery status. Whenever this is low, your stimulator should be charged as soon as possible.
- **Keep your base station plugged in** so it is **always ready to charge** your charger when needed.



This document is not intended as a substitute for Spectra WaveWriter™ System Patient Manuals. The contents of the handbook must be carefully read and fully understood before using your Charging System. If you have any questions, or need clarification on anything contained in this manual, please [contact your physician](#).

CAUTION: Indications, contraindications, warnings and instructions for use can be found in the product labelling supplied with each device.

Information for the use only in countries with applicable health authority product registrations.

Illustrations for information purposes – not indicative of actual size or clinical outcome.

NM-608106-AA DINNM0160EA



Boston Scientific Pty Ltd

PO Box 332

Botany NSW 1455

Australia

Tel +61 2 8063 8100

Fax +61 2 9330 1404

© 2019 Boston Scientific Corporation
or its affiliates. All rights reserved.

ANZ_PSST_19233_AA AUG 2019